

LiNA Versascope™ Sheath

Single-use Sheath for Diagnostic and Operative Hysteroscopy



Enables See & Treat Hysteroscopy



Innovation in Gynecology

LiNA Versascope™ Sheath

Single-use Sheath for Diagnostic and Operative Hysteroscopy

See and Treat Hysteroscopy

Minimally invasive for patient comfort

- Small diameter ensures ease of insertion, often without dilation.
- Minimizes the need for general anesthesia.

Effective distension and clear visualisation

- The LiNA Versascope™ Sheath features inflow and outflow channels, for optimal fluid management.


Versatility

- The expandable working channel allows the use of a wide range of instruments.




Single-use Instrumentation


- Our affordable selection of single-use instrumentation enables see and treat hysteroscopy without the limitations of traditional reusable instruments.

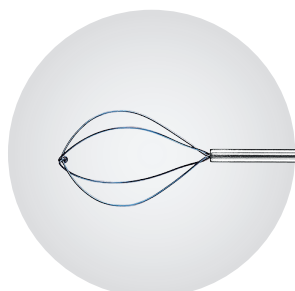


Lasso

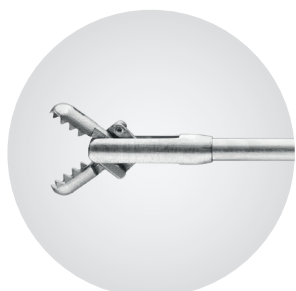

Biopsy Forceps

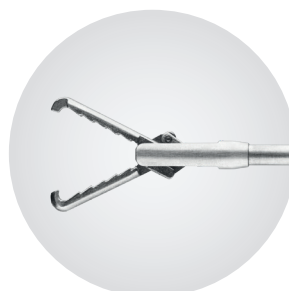


Rat Tooth
Alligator Grasper

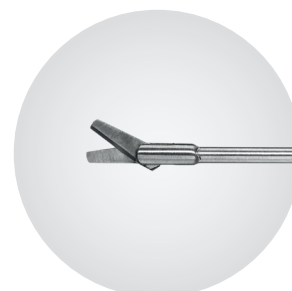

Scissors

Fully Integrated for Convenience & Utility

The expandable LiNA Versascope™ Sheath is compatible with the 1.9 mm AlphaScope™ to form an integrated system designed to enable quick diagnostic-to-therapeutic conversions during hysteroscopic procedures.

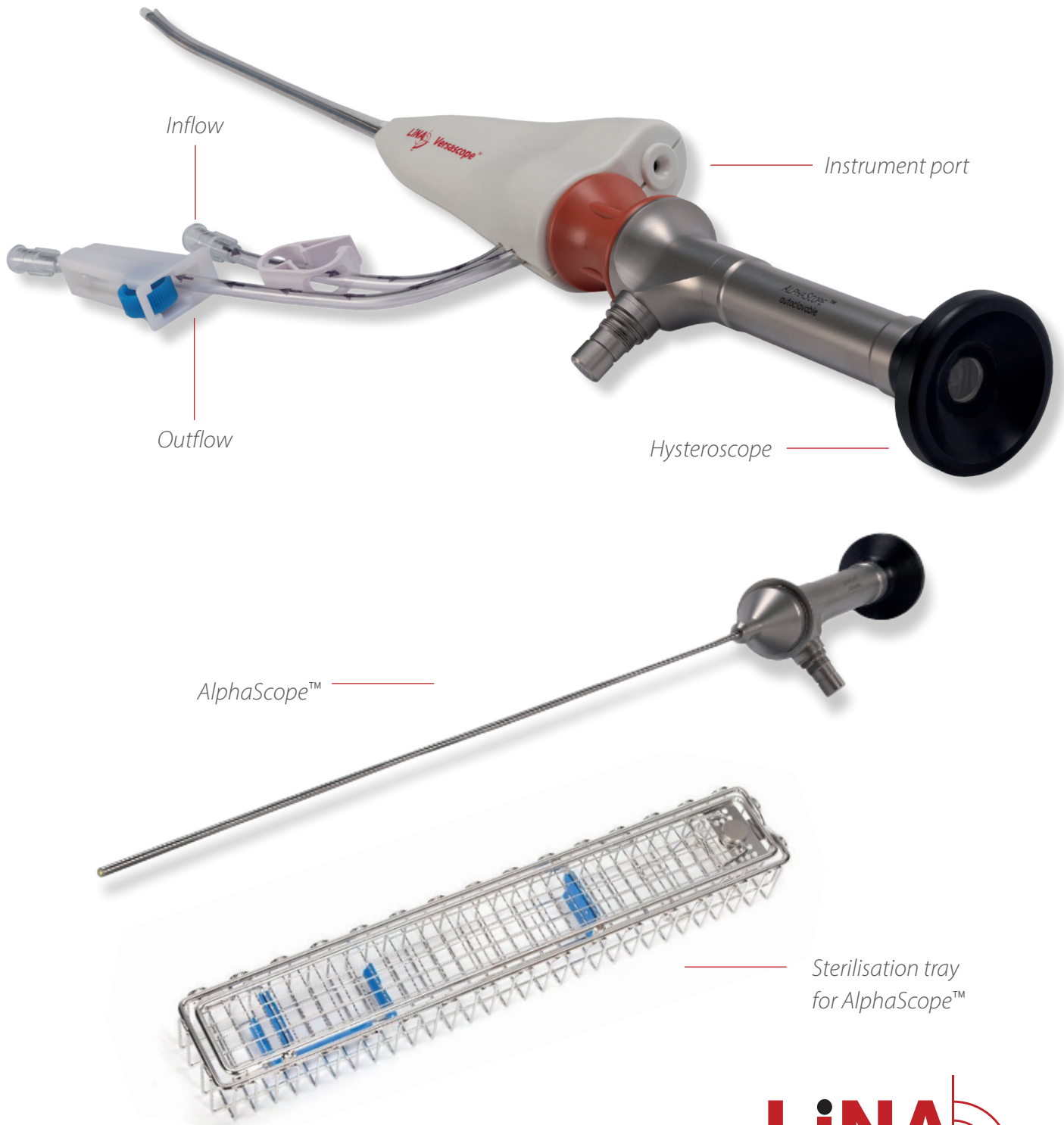

LiNA Versascope™ Sheath

Single-use Sheath for Diagnostic and Operative Hysteroscopy

Technical Specifications:

LiNA Versascope™ Sheath	
Maximum Diameter (diagnostic)	3.7 mm
Expandable Working Channel	2 mm
Instrument channel working length	280 mm

Gimmi® AlphaScope™	
Working length/ Scope length	280 mm / 410,4 mm
Diameter	1,85 mm
Directions of view / Field of view	$\alpha=0^\circ / \beta=80^\circ$
Pixel	40.000 Pixel
Ocular	0° straight view
Light conductor connection	fixed: ACMI; additional: Storz/Olympus and Wolf
Material	Scope body + sheath, Light post: stainless steel Eyepiece housing: PEEK Distal and proximal; window: sapphire glass Distal and proximal glue: Expoxide, Polyepoxide

Ordering Information:

 Gimmi® AlphaScope™ 0123

Product code	Description	UOM
GMS40A	AlphaScope™ 1.9 mm Hysteroscope	Box of 1ea
GMS710A	Tray for AlphaScope™ Hysteroscope	Box of 1ea
GMS773A	2,3 m Fiber Optic Light Cable including 4 Adapter set	Box of 1ea
GMS774A	3 m Fiber Optic Light Cable including 4 Adapter set	Box of 1ea

Ordering Information:

 LiNA Versascope™ 0197

Product code	Description	UOM
GMS805	LiNA Versascope™ Sheath	Box of 5ea
GMS720	LiNA Versascope™ Outflow Cannula	Box of 10ea

Ordering Information:

 LiNA Accessories 2274

Product code	Description	UOM
OP-LAS10-3-I	LiNA OperåScope™ Lasso (10 mm, sterile single-use)	Box of 3ea
OP-LAS16-3-I	LiNA OperåScope™ Lasso (16 mm, sterile single-use)	Box of 3ea
OP-FCP-6-I	LiNA OperåScope™ Biopsy Forceps (sterile single-use)	Box of 6ea
OP-GRA-6-I	LiNA OperåScope™ Rat Tooth Alligator Grasper (sterile single-use)	Box of 6ea
OP-SCI-3-I	LiNA OperåScope™ Scissors (sterile single-use)	Box of 3ea
OP-TDK-6-I	LiNA OperåScope™ Tubing & Drape Kit (sterile single-use)	Box of 6ea

©2026 LiNA Medical AG. All rights reserved. All products may not be available for sale in all geographies. Please contact your local LiNA Medical representative for information about products available in your area. This information is intended for healthcare professionals. If you are a patient seeking more information, please consult your healthcare provider. Before using any LiNA product, always refer to the IFU included in our packaging for complete indications for use, contraindications, warnings, precautions, and instructions. LiNA is a registered trademark of LiNA Medical in EU, US and CN.

EU MDR 2017/745 device class for REF:

- GMS805 and GMS720: class Is

Name and number of the Notified Body: TÜV Rheinland LGA Products GmbH #0197

- OP-SCI-3-I, OP-GRA-6-I, OP-LAS10-3-I, OP-LAS16-3-I, OP-FCP-6-I: class Is

Name and number of the Notified Body: TÜV Nord Polska Sp. z o.o. # 2274.

- GMS40A: class IIa

- GMS710A, GMS773A: class I

Name and number of the Notified Body: TÜV SUD Product Service GmbH #0123



LiNA Medical AG, Platz 3, 6039 Root D4 | Switzerland

Phone: +41 78 617 56 56 , info@lina-medical.com | www.linamed.com



LiNA Medical Polska Sp. z o.o., Rolna 8A, Sady, 62-080 Tarnowo Podgórze | Poland

Phone: +48 61 222 21 21, info.poland@lina-medical.com | www.lina-medical.pl

EN ISO 13485:2016

Quality Management
System

